

IoT-Line Precision Balance KERN 572



- Protective working cover included with delivery

Technical data

- Large backlit LCD display, Digit height 21 mm
- Dimensions weighing surface, stainless steel
 - A** \varnothing 106 mm
 - B** \varnothing 150 mm
 - C** WxD 160x200 mm, see larger picture
- Overall dimensions WxDxH 180x310x85 mm
- Net weight **A**, **B** approx. 2,4 kg **C** approx. 2,8 kg
- Permissible ambient temperature -10 °C/40 °C


Accessories

- Protective working cover, scope of delivery: 5 items, KERN 572-A02S05
- Internal rechargeable battery pack, operating time up to 48 h without backlight, charging time approx. 8 h, KERN YKR-01
- External data interface RS-232, interface cable included, KERN KUP-01
- External data interface USB, interface cable included, KERN KUP-03
- Interface adapter Ethernet, KERN KUP-04
- Interface adapter WiFi, KERN KUP-05
- Bluetooth interface adapter, KERN KUP-06
- Extension box for connecting up to three interfaces in parallel, KERN KUP-13
- **B**, **C** Hook for underfloor weighing, KERN 572-A03
- **A** Large glass draught shield with 3 sliding doors for easy access to the items being weighed. Weighing space WxDxH 150x140x130 mm, for models with weighing plate size **A**, KERN 572-A05

All-rounder, e.g. as precision balance in the laboratory or in harsh industrial applications, ideal for the diverse possibilities of Industry 4.0 applications

Features

- Thanks to the many typical laboratory functions, such as, for example, recipe function, percentage determination, combined with the high level of precision, the KERN 572 is a reliable partner for day-to-day work in the laboratory
- The robust version, typical industrial functions, such as piece-counting, vibrationfree weighing and the large weighing ranges also make these balances ideal for all industrial applications, where a high level of precision is required
- KERN Universal Port (KUP): permits the connection of an external KUP interface adapter, such as, for example, RS-232, USB, Bluetooth, WiFi or Ethernet, for the exchange of data and control commands, without any installation outlay
- KERN Communication Protocol (KCP): The KCP permits searching and remote control of the balance using external control devices or computers
- For further information on KUP and KCP see page 20/21
- Standardised, simplified concept of operation
- Freely programmable weighing unit, e.g. display direct in special units such as length of thread g/m, paper weight g/m², or similar
- The robust aluminium diecast housing maintains the stability, protects the weighing technology elements and is robust enough to cope with everyday use
- **A** Ring-shaped draught shield standard, only for models with weighing plate size **A**, weighing space \varnothing xH 157x43 mm
- Level indicator to level the balance precisely
- Underfloor weighing: load support on the underside of the balance for models. Loop for underfloor weighing, standard for models with weighing plate size **A**

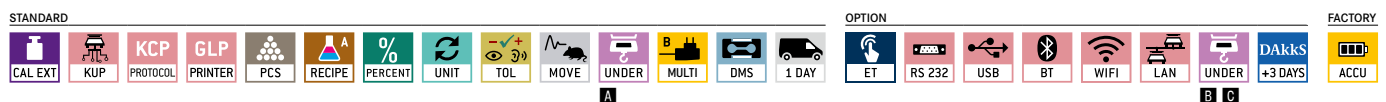

Model	Weighing capacity [Max]	Readability	Reproducibility	Linearity	Resolution	Weighing plate	Options
KERN	g	[d] g	g	g	Points		DAkkS Calibr. Certificate DAkkS KERN
572-30	240	0,001	0,001	± 0,003	240.000	A	963-127
572-31	300	0,001	0,002	± 0,005	300.000	A	963-127
572-32	420	0,001	0,002	± 0,005	420.000	A	963-127
572-33	1600	0,01	0,01	± 0,03	160.000	B	963-127
572-35	2400	0,01	0,01	± 0,03	240.000	B	963-127
572-37	3000	0,01	0,02	± 0,05	300.000	B	963-127
572-39	4200	0,01	0,02	± 0,05	420.000	B	963-127
572-45	12000	0,05	0,05	± 0,15	240.000	C	963-128
572-55	20000	0,05	0,1	± 0,25	400.000	C	963-128
572-43	10000	0,1	0,1	± 0,3	100.000	C	963-128
572-49	16000	0,1	0,1	± 0,3	160.000	C	963-128
572-57	24000	0,1	0,1	± 0,3	240.000	C	963-128

CAL INT
Internal adjusting
 Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)

CAL EXT
Adjusting program CAL
 For quick setting up of the balance's accuracy. External adjusting weight required

ET
EasyTouch
 Suitable for the connection, data transmission and control through PC or tablet

MEMORY
Memory
 Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.

ALIBI
Alibi memory
 Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.

KUP
KERN Universal Port (KUP)
 allows the connection of external KUP interface adapters, e.g. RS-232, RS-485, SB, Bluetooth, WIFI, Analogue, Ethernet etc. for the exchange of data and control commands, without installation effort

RS 232
RS-232 Data interface
 To connect the balance to a printer, PC or network

RS 485
RS-485 Data interface
 To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible

USB
USB Data interface
 To connect the balance to a printer, PC or other peripherals

BT
Bluetooth* Data interface
 To transfer data from the balance to a printer, PC or other peripherals

WIFI
WIFI Data interface
 To transfer data from the balance to a printer, PC or other peripherals

SWITCH
Control outputs
 (optocoupler, digital I/O)
 To connect relays, signal lamps, valves, etc.

ANALOG
Analogue interface
 to connect a suitable peripheral device for analogue processing of the measurements

DUAL
Interface for second balance
 For direct connection of a second balance

LAN
Network interface
 For connecting the scale to an Ethernet network

KCP PROTOCOL
KERN Communication Protocol (KCP)
 It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems

GLP INTERN
GLP/ISO log intern
 The balance displays weight, date and time, independent of a printer connection

GLP PRINTER
GLP/ISO log Printer
 With weight, date and time. Only with KERN printers.

PCS
Piece counting
 Reference quantities selectable. Display can be switched from piece to weight

RECIPE A
Recipe level A
 The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out

RECIPE B
Recipe level B
 Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display

SUM A
Totalising level A
 The weights of similar items can be added together and the total can be printed out

PERCENT
Percentage determination
 Determining the deviation in % from the target value (100 %)

UNIT
Weighing units
 Can be switched to e.g. nonmetric units. See balance model. Please refer to KERN's website for more details

TOL
Weighing with tolerance range (Checkweighing)
 Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model

MOVE
Hold function
 (Animal weighing program)
 When the weighing conditions are unstable, a stable weight is calculated as an average value

IP
Protection against dust and water splashes IPxx
 The type of protection is shown in the pictogram

UNDER
Suspended weighing
 Load support with hook on the underside of the balance

BATT
Battery operation
 Ready for battery operation. The battery type is specified for each device

ACCU
Rechargeable battery pack
 Rechargeable set

MULTI
Universal plug-in power supply
 with universal input and optional input socket adapters for
 A) EU, CH, GB
 B) EU, CH, GB, US
 C) EU, CH, GB, US, AUS

230 V
Plug-in power supply
 230V/50Hz in standard version for EU, CH. On request GB, USA or AUS version available

230 V
Integrated power supply unit
 Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request

DMS
Weighing principle Strain gauges
 Electrical resistor on an elastic deforming body

T-FORK
Weighing principle Tuning fork
 A resonating body is electromagnetically excited, causing it to oscillate

FORCE
Weighing principle Electromagnetic force compensation
 Coil inside a permanent magnet. For the most accurate weighings

SC TECH
Weighing principle Single cell technology
 Advanced version of the force compensation principle with the highest level of precision

M +3 DAYS
Conformity Assessment
 The time required for conformity assessment is specified in the pictogram

DAkkS +3 DAYS
DAkkS calibration possible (DKD)
 The time required for DAkkS calibration is shown in days in the pictogram

ISO +4 DAYS
Factory calibration (ISO)
 The time required for Factory calibration is shown in days in the pictogram

1 DAY
Package shipment
 The time required for internal shipping preparations is shown in days in the pictogram

2 DAYS
Pallet shipment
 The time required for internal shipping preparations is shown in days in the pictogram

* The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.